Bilancia Parlante

Talking scale - Báscula hablante - Pèse-personne vocal

ISTRUZIONI

INSTRUCTIONS

INSTRUCCIONES

INSTRUCTIONS

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INSTRUCTIONS

USE AND ACTIVATION

- Insert the supplied batteries (with the correct polarity) paying attention to the + and - signs in the battery compartment. If the product does not light up, remove the batteries and reinsert them.
- For the first use (or after each movement), always calibrate the scale for maximum precision using the procedure indicated in the quick guide inside the package.
- Always place the scale on a hard, level surface to ensure accuracy.

HOW TO CHANGE THE LANGUAGE

To change the language, press the **COMMAND** button located behind the scale only once, choosing one off the four available.

ITALIAN

ENGLISH

SPANISH

FRENCH

HOW TO SILENCE OR ACTIVATE THE VOICE

To silence or activate the voice, press the **COMMAND** button behind the scale for 3 seconds.

ACTIVATED VOICE

ITEM DEACTIVATED


OVERLOAD INDICATOR

When the scale is overloaded (over 180 kg), the following message will be displayed on the screen:

BATTERY REPLACEMENT INDICATOR

When the batteries need to be replaced, the following message will appear on the screen:

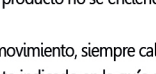CONTENTS OF THE PACKAGE

- 1 Weighing scale
- 2 1.5 V AAA batteries
- 1 Instruction manual
- 1 Quick Guide

TECHNICAL DATA

Model: iVoice
Weighing surface size 30 x 30 cm
Maximum capacity: 180 Kg
Weight subdivision: 0.1 Kg
1.5 V batteries type AAA
Material: Glass, plastic, metal

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The product is covered by a 2-year warranty starting from the date of purchase, in case of defects of the material and manufacturing of the scales. The guarantee does not apply:

- In case of damage caused by improper use.
- For components subject to wear.
- In case of responsibility of the customer.

- If the appliance was opened at an unauthorized workshop.

The rights of the purchaser under the law arising from the provision of the guarantee are not affected by this warranty.

To claim the warranty within its validity period, the customer must prove the purchase of the product. It is possible to assert your warranty claims against VIBAS within a period of 2 years.

In case of complaints, contact our service department at the following contacts:

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